

AML-193 Human Peripheral Blood Acute Monocyte Leukemia Cell Line

SKU: CDS4063 | Category: Cell Lines — Cancer Cell Lines

Validation & Biomarker Data Guide · Lineage classification: Monocyte / macrophage

1. Product Overview

AML-193 Human Peripheral Blood Acute Monocyte Leukemia Cell Line is an established tumor-derived cell line for in vitro oncology research. It is suitable for studies of tumor cell biology, proliferation and viability, signaling, biomarkers, drug response, cytotoxicity and cell-based assay development. The model can be incorporated into screening and mechanistic workflows under validated culture conditions. For research use only.

Growth type: Suspension Cells | Tumorigenic status: Tumor Cells | Pack size: 1 x 1 million cells

2. Typical Applications

- Cell-based assay development, pathway and mechanism studies, biomarker characterization, and comparative functional testing consistent with this lineage.

3. Identity & Phenotypic Marker Panel

A representative identity panel for monocyte / macrophage products is shown below. Actual per-lot Certificates of Analysis should be consulted for verified identity data.

Marker	Expected Result	Range
CD14	Positive	75–97% positive
CD68	Positive	70–95% positive
CD11b	Positive	70–95% positive
CD3	Negative / low	≤5% (background)
CD19	Negative / low	≤4% (background)

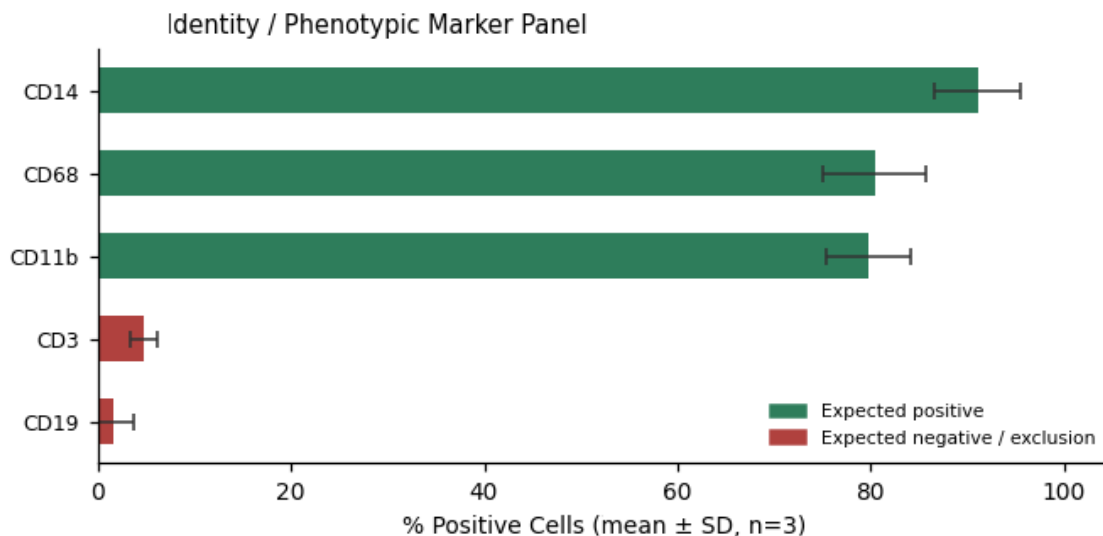


Figure 1. Marker expression profile (mean $\pm$ SD, n = 3 replicates).

4. Growth Kinetics & Viability

Growth and viability trends over 120 hours in standard culture conditions.

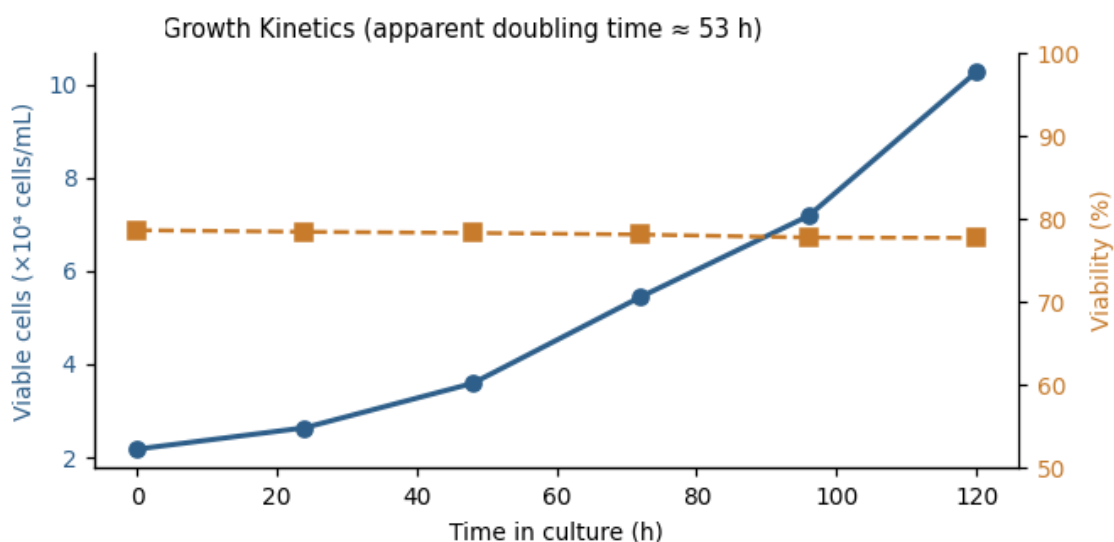


Figure 2. Growth curve (viable cell density, left axis) and viability trend (right axis).

5. Functional Validation Assay

Phagocytosis assay: Uptake of fluorescently labeled latex beads or bioparticles quantified by flow cytometry/imaging to confirm phagocytic competence; LPS-stimulated TNF- α /IL-6 secretion by ELISA.

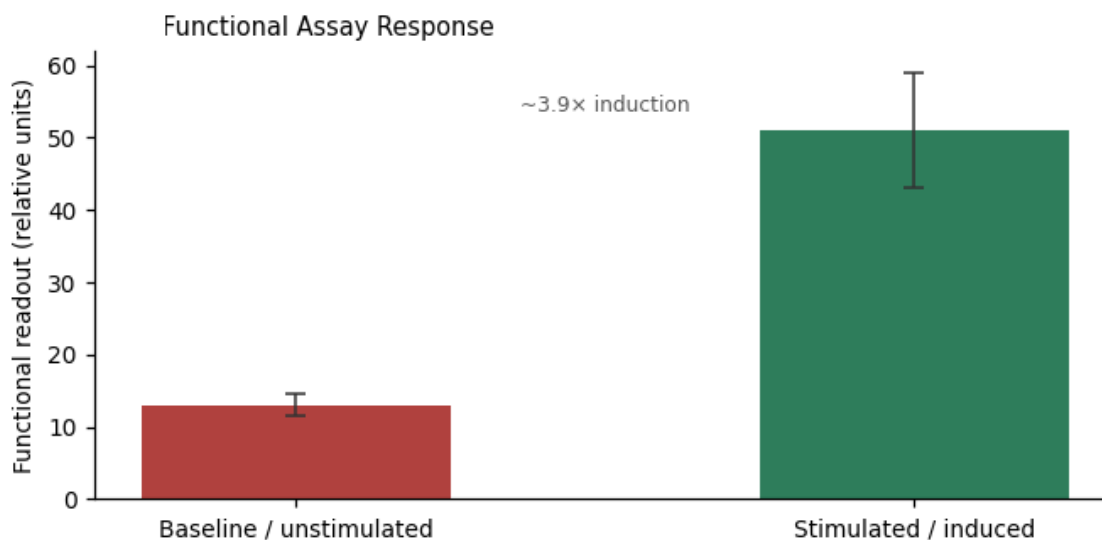


Figure 3. Functional assay readout — ~3.9× induction over baseline.

6. QC Acceptance Criteria

Acceptance Criterion	Result
≥70% CD14+/CD68+ by flow cytometry or IHC	PASS
Positive phagocytic uptake (>60% bead-positive cells)	PASS
Apparent population doubling time	≈ 53 hours
Functional assay outcome	~3.9× induction over baseline

7. Handling & Culture Notes

- Refer to the product-specific Certificate of Analysis and datasheet for validated culture medium, seeding density, and passaging recommendations.
- As with any primary cell or cell line, confirm mycoplasma-negative status, expected morphology by phase-contrast microscopy, and viability by trypan blue or equivalent dye-exclusion assay at each passage before use in downstream assays.