

# HMEC-1 Human Microvascular Endothelial Cell Line

SKU: CDS4852 | Category: Cell Lines — Normal and Non-Tumor Cell Lines

**Validation & Biomarker Data Guide** · Lineage classification: Endothelial cell

## 1. Product Overview

HMEC-1 Human Microvascular Endothelial Cell Line is an established non-tumor cell model for in vitro research. It is suitable for studies of cell physiology, growth, signaling, toxicity, drug response, biomarker analysis and assay development under controlled culture conditions. The cell line provides a reproducible system for basic and translational cell biology. For research use only.

Growth type: Adherent Cells | Tumorigenic status: Normal Cells | Pack size: 1 x 1 million cells

## 2. Typical Applications

- Cell-based assay development, pathway and mechanism studies, biomarker characterization, and comparative functional testing consistent with this lineage.

## 3. Identity & Phenotypic Marker Panel

A representative identity panel for endothelial cell products is shown below. Actual per-lot Certificates of Analysis should be consulted for verified identity data.

| Marker          | Expected Result | Range             |
|-----------------|-----------------|-------------------|
| CD31 (PECAM-1)  | Positive        | 85–99% positive   |
| VE-cadherin     | Positive        | 75–97% positive   |
| vWF             | Positive        | 70–95% positive   |
| Pan-cytokeratin | Negative / low  | ≤5% (background)  |
| α-SMA           | Negative / low  | ≤10% (background) |

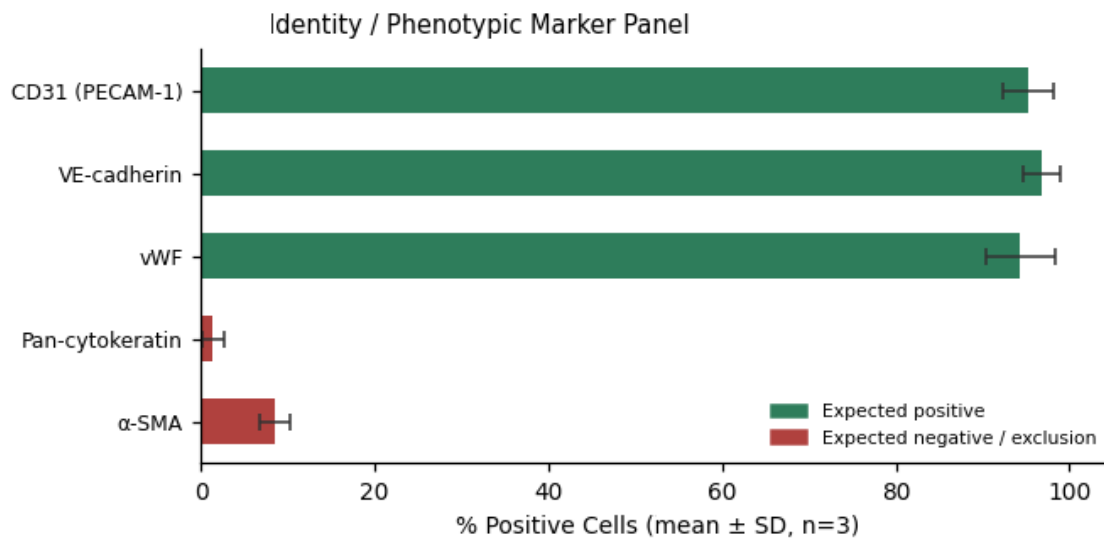


Figure 1. Marker expression profile (mean ± SD, n = 3 replicates).

#### 4. Growth Kinetics & Viability

Growth and viability trends over 120 hours in standard culture conditions.

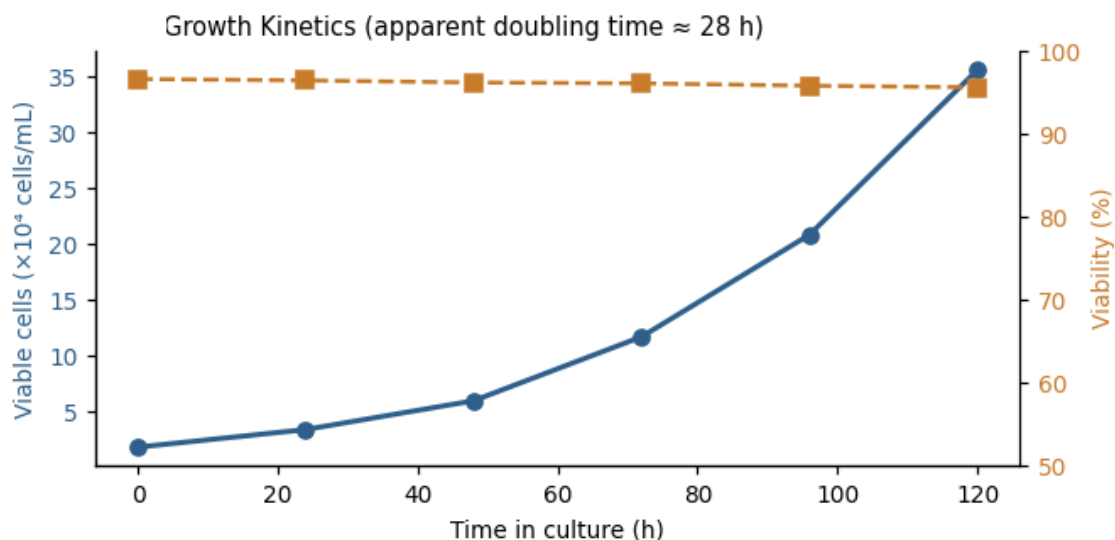


Figure 2. Growth curve (viable cell density, left axis) and viability trend (right axis).

#### 5. Functional Validation Assay

**Tube formation (angiogenesis) assay:** Matrigel tube-formation assay quantifying network/branch-point formation to confirm angiogenic functional competence; Ac-LDL uptake as secondary marker.

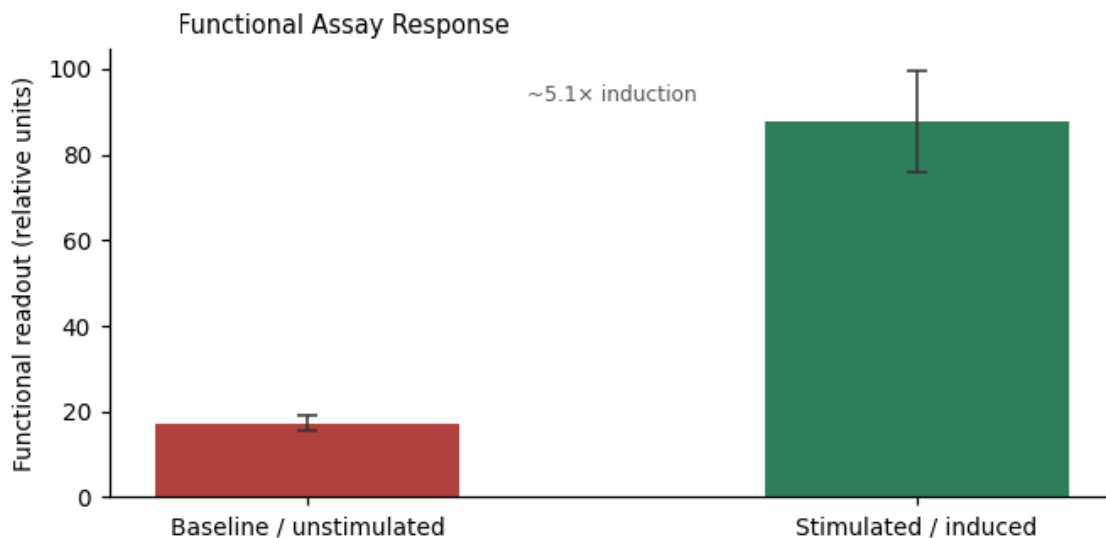


Figure 3. Functional assay readout — ~5.1× induction over baseline.

## 6. QC Acceptance Criteria

| Acceptance Criterion                              | Result                        |
|---------------------------------------------------|-------------------------------|
| ≥85% CD31+ by flow cytometry/ICC                  | PASS                          |
| Positive tube formation on Matrigel within 6-18 h | PASS                          |
| Apparent population doubling time                 | ≈ 28 hours                    |
| Functional assay outcome                          | ~5.1× induction over baseline |

## 7. Handling & Culture Notes

- Refer to the product-specific Certificate of Analysis and datasheet for validated culture medium, seeding density, and passaging recommendations.
- As with any primary cell or cell line, confirm mycoplasma-negative status, expected morphology by phase-contrast microscopy, and viability by trypan blue or equivalent dye-exclusion assay at each passage before use in downstream assays.