

## Immortalized Cell Line Of The Human Ureteric Epithelium, #SV-HUC-1


Immortalized Cell Line Of The Human Ureteric Epithelium, #SV-HUC-1 is an immortalized cell model for sustained in vitro investigation of the named cell type.

Cell-based models help researchers examine biological processes under controlled laboratory conditions.

They offer a practical way to relate environmental or treatment changes to cellular growth and function.

|                     |                                 |
|---------------------|---------------------------------|
| CTDS SKU            | CDS10199                        |
| Pack size           | 1 x 1 million cells             |
| Product category    | Cell Lines                      |
| Product subcategory | Normal and Non-Tumor Cell Lines |

### PRODUCT IMAGE

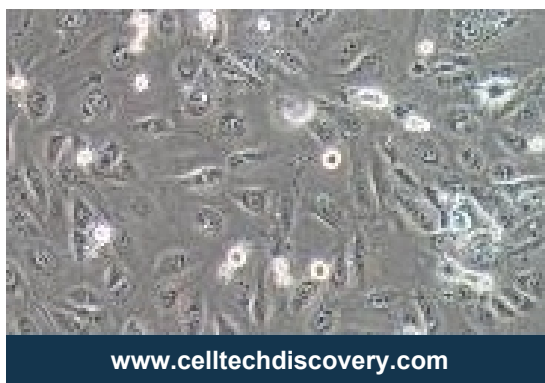

### PRODUCT BACKGROUND

Immortalized Cell Line Of The Human Ureteric Epithelium, #SV-HUC-1 is an established non-tumor cell model for in vitro research. It is suitable for studies of cell physiology, growth, signaling, toxicity, drug response, biomarker analysis and assay development under controlled culture conditions. The cell line provides a reproducible system for basic and translational cell biology. For research use only.

### APPLICATIONS

#### 1. Cancer Biology Research

Mechanistic Studies

Gene Function Studies

Signal Transduction

#### 2. Drug Discovery and Development

Screening of Anticancer Compounds

IC50/EC50 Studies

Mechanism of Action

Combination Therapies

#### 3. Assay Development

ELISA and Immunoassays

Flow Cytometry Assays

Reporter Gene Assays

#### 4. Omics and Biomarker Research

Transcriptomics / Genomics

Proteomics

Metabolomics

## 5. Tumor Heterogeneity and Resistance Mechanisms

Clonal Selection

Acquired Resistance Models

## 6. Personalized Medicine and Precision Oncology

Patient-Derived Cell Lines

Companion Diagnostics

## 7. Immuno-oncology

Co-culture with Immune Cells

Checkpoint Inhibitor Testing

Antibody-Dependent Cytotoxicity (ADCC)

## 8. Vaccine Development

Neoantigen Discovery

Antigen Presentation Studies

## BIOLOGICAL CHARACTERISTICS

|                              |                           |
|------------------------------|---------------------------|
| Alternate name 1             | CVCL_1786 Cells           |
| Alternate name 2             | sv-huc1 Cells             |
| Tissue / origin              | Ureteric Epithelium       |
| Cell growth type             | Mixed Adherent/Suspension |
| Biosafety level              | BSL-2                     |
| Tumorigenicity / cell status | Normal Cells              |

## CULTURE REQUIREMENTS

|                 |                                                 |
|-----------------|-------------------------------------------------|
| Culture medium  | F12K + 10% FBS + 1%PS                           |
| Freezing medium | 55% Complete Culture Medium + 40% FBS + 5% DMSO |

## RECEIPT AND CRYOPRESERVATION

Upon Receipt of Cells

1 mL cryovials are packed in dry ice for transportation. After receipt, clean the cryovial with 70% IPA to disinfect. Stored in a -80°C refrigerator overnight and then transfer to liquid nitrogen or revive them directly as per your use.

Note If the dry ice is found to have evaporated or the cryovial cap is detached/ damaged or the cells seems damaged, please contact us immediately.

## CULTURE AND RECOVERY INSTRUCTIONS

Recovery of Frozen Cells:

For Suspension Cell passaging, the following methods can be used as reference:

- 1) Shake the cryotube containing 1mL of cell suspension quickly within seconds in a 37°C water bath to thaw.
- 2) Add it to a centrifuge tube containing 4-6mL of complete medium and mix well. Centrifuge at 1000rpm for 3-5mins.
- 3) Discard the supernatant, and resuspend the cells in complete medium.
- 4) Then add the cell suspension to a culture flask (or dish) containing 6-8ml of complete medium and culture overnight at 37°C.
- 5) Observe the cell growth and cell density under a microscope the next day.
- 6) Cell subculture: If the cell density reaches 80%-90%, subculture can be performed.

For Adherent Cell passaging, the following methods can be used as reference:

- 1) Shake the cryotube containing 1mL of cell suspension quickly within seconds in a 37°C water bath to thaw.
- 2) Add it to a centrifuge tube containing 4-6mL of complete medium and mix well. Centrifuge at 1000rpm for 3-5mins.
- 3) Discard the supernatant, and rinse the cells 1-2 times with PBS without calcium and magnesium ions.
- 4) Add 0.25% (w/v) trypsin-0.53 mM EDTA to the culture flask (1-2 mL for T25 flask, 2-3 mL for T75 flask),
- 5) Place in a 37°C incubator for digestion for 1-2 minutes (difficult-to-digest cells can be properly prolonged for digestion time),
- 6) Observe the cell digestion under a microscope. If most of the cells become round and fall off, quickly take it back to the

operating table, tap the culture flask a few times, and add 3-4 ml of culture medium containing 10% FBS to terminate digestion.

7) Gently mix and aspirate, centrifuge at 1000rpm for 3-5min, discard the supernatant. Add 1-2mL culture medium and mix evenly.

8) Divide the cell suspension into new T25 bottles at a ratio of 1:2.

9) Add 6-8ml of new complete culture medium prepared according to the instructions to maintain cell growth vitality, and perform subsequent passages at a ratio of 1:2~1:5 as per actual conditions.

## CHARACTERIZATION AND QUALITY INFORMATION

|                              |                                                                                                                                                                                                                                |
|------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Sterility                    | Bacteria: Negative                                                                                                                                                                                                             |
|                              | Yeast: Negative                                                                                                                                                                                                                |
|                              | Mycoplasma: Negative                                                                                                                                                                                                           |
| Pathogen testing             | HIV: Negative                                                                                                                                                                                                                  |
|                              | Hepatitis B: Negative                                                                                                                                                                                                          |
| Markers / assay readout      | Phenotype QC: viability dye plus 2-3 lineage-positive markers and at least one exclusion marker. For human continuous lines, use STR profiling for identity; flow Phenotype is supportive only.                                |
| Assay / validation method    | Flow cytometry marker characterization                                                                                                                                                                                         |
| Authentication / QC evidence | Phenotype QC: viability dye plus 2-3 lineage-positive markers and at least one exclusion marker. For human continuous lines, use STR profiling for identity; flow Phenotype is supportive only.; Sterility: Bacteria: Negative |
|                              | Yeast: Negative                                                                                                                                                                                                                |
|                              | Mycoplasma: Negative; Pathogen testing: HIV: Negative Hepatitis B: Negative; Parental-line authentication: confirm applicability to this product                                                                               |

## HANDLING AND STORAGE

|              |                                                          |
|--------------|----------------------------------------------------------|
| Shelf life   | 12 months at time of shipping                            |
| Shipping     | Shipped on Dry Ice at -20°C                              |
| Storage      | Store at -176°C.                                         |
| Intended use | For Research Use Only. Not for diagnostics or human use. |

## PRODUCT DISCLAIMER

The data indicated herein are as indicated and validated in our laboratory. As these cell lines are of human origin, they carry a risk of latent human viruses (e.g., Epstein-Barr virus, hepatitis B). Therefore they should be handled accordingly with precaution. These reagents are for research use only and not for in-vitro diagnostics or human use.

## ADDITIONAL CLASSIFICATION

|                  |                               |
|------------------|-------------------------------|
| Product grouping | Normal / Non-Tumor Cell Lines |
| Product subgroup | Human                         |

Not a lot-specific certificate of analysis.